



# ConMed Suction Coagulators


- Improved ergonomics
- Malleable cannula for easy customization
- Both handswitch and footswitch models available



# Suction Coagulation

Electrosurgery has been proven essential for surgeons in controlling bleeding. Suction coagulators provide a means of aspirating blood and small tissue particulate in addition to controlling bleeding during open surgical procedures.

With improved ergonomics and a malleable cannula, ConMed suction coagulators are easily customizable to suit any surgeon's needs.

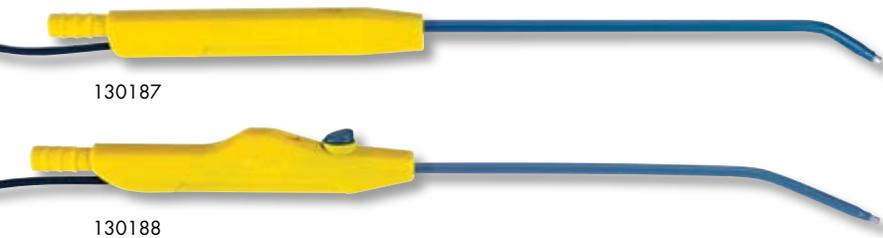

## ConMed Suction Coagulators Technical Information

Type	Footswitch or handswitch
Suction application	Place finger over suction vent on handle
Handle construction	Ergonomically designed yellow ABS
Cannula	Aluminum alloy
Cannula coating	Blue polyolefin
Obturator	Polypropylene cleaning tube and protective insert for bending
Suction tubing	0.6cm (1/4 in.) or 0.5cm (3/16 in.) tubing (not included) compatible
Cordset	Footswitching– 10' (3.0m) Handswitching – 10' (3.0m)

## Ordering Information

Catalog Number	Product Description	Packaging
130187	Suction Coagulator, Foot Control, 10 French, 6" (15 cm), Monopolar, Single Use; Sterile	25/case 1/pkg
130188	Suction Coagulator, Hand Control, 10 French, 6" (15 cm), Monopolar, Single Use; Sterile	25/case 1/pkg

Please contact your local ConMed Electrosurgery Representative or Customer Service at 1-800-448-6506 today to arrange an evaluation.

### CONMED CORPORATION PRODUCT AREAS:

ARTHROSCOPY • ELECTROSURGERY • ENDOSCOPY • ENDOSURGERY • GASTROENTEROLOGY • INTEGRATED SYSTEMS • PATIENT CARE • POWERED INSTRUMENTS • PULMONOLOGY



525 French Road  
 Utica, NY 13502  
 (315) 797-8375  
 Customer Service: 1-800-448-6506  
 FAX: 1-800-438-3051  
 International Sales: +1(315) 797-8375 • FAX: +1(315) 735-6235  
 email: info@mail.conmed.com  
 www.conmed.com