



solace
IN-FLOOR DIALYSIS SCALES




IN-FLOOR DIALYSIS PATIENT WEIGHING

- Flush-Mounted with Floor for Easy Accessibility
- Stainless Steel, Anti-Corrosion Weighing Platform
- Multi-Color, Touchscreen Indicator/Printer Combo Included



DETECTO's USA-made Solace series are available in 6 x 4 ft sizes for weighing stretchers.

DETECTO's Solace series flush-mounted platform scales eliminate tripping hazards and hallway accessibility issues.

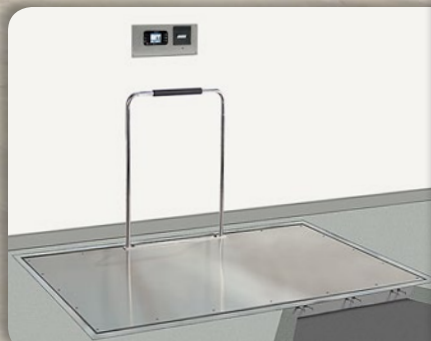



IN-FLOOR DIALYSIS SCALES

IN-FLOOR DIALYSIS SCALES

DETECTO's Solace series USA-made in-floor dialysis scales are flush-mounted platform scales that feature stainless steel decks and include a pit frame and trim ring. The flush surface eliminates tripping hazards and hallway accessibility issues.

The Solace series is available with or without a handrail and all models include a recessed-mount indicator and printer combination. These in-floor scales easily accommodate patients with walkers, canes, wheelchairs, or seated in dialysis chairs.



The ID models are designed to be mounted in a pit flush with the floor. A tubular handrail with padded handgrip is available on some models for walk-on patients.

DIALYSIS + BARIATRIC WEIGHING

- 1,000 lb x 0.2 lb / 450 kg x 0.1 kg capacity
- 8 standard models to choose from
- Handrails available for patient stability while weighing
- Indicator/printer, pit frame, scale, and trim ring included
- Standard platform sizes: 3 x 3 ft, 4 x 3 ft, 4 x 4 ft, and 6 x 4 ft
- Custom sizes and configurations are available
- USA-made quality construction



RECESSED-MOUNT INDICATOR + PRINTER COMBO INCLUDED



The integrated scale displays on DETECTO's USA-made multi-color indicator that is wall-mounted with a stainless steel cover making it easy for clinicians to read and accurately record measurements.

DETECTO's Solace series in-floor scales are covered with durable stainless steel to provide years of maintenance-free use. The built-in thermal printer makes it easy to print out a weight record for patient recordkeeping.

- Multi-color, touchscreen LCD display
- Built-in thermal ticket printer (included)
- Displays weight, height, and BMI onscreen simultaneously
- Large 0.9-inch-high, clinical-blue LCD weight digits
- Push button or manual entry of tare weights for wheelchairs and stretchers
- Automatically locks on weight temporarily for a few seconds in order to easily view and record measurement
- Menu button provides quick access to user-friendly menu settings for scale customization
- Standard 110V AC power 110V – 240V AC universal switching power supply

ID: 5432
GROSS: 489.0 lb
TARE: 100.0 lb
NET: 389.0 lb

ID: 8976
GROSS: 218.2 lb
TARE: 50.0 lb
NET: 168.2 lb

The Solace series in-floor scale's ticket printer can print Gross, Tare, Net and Patient IDs for record keeping.



The Solace series is flush mounted with the floor and includes the pit frame and trim ring, plus an interface cable from the scale to indicator.



DETECTO's Solace series in-floor scales are ideal for dialysis clinics, emergency rooms, and clinical facilities where fast, efficient patient measurement is necessary.

Complete in-floor scale package for continuous patient weight monitoring



In-floor scales comes with model 855RMP digital weight indicator/printer combo, pit frame, scale, stainless steel deck, and trim ring. Custom sizes and configurations are available.



SOLACE SERIES IN-FLOOR DIALYSIS SCALES



Model	ID-3636S-855RMP	ID-3636SH-855RMP	ID-4836S-855RMP	ID-4836SH-855RMP
Capacity	1,000 lb x .2 lb / 450 kg x .1 kg			
Platform Size	36 x 36 in / 91 x 91 cm	36 x 36 in / 91 x 91 cm	48 x 36 in / 122 x 91 cm	48 x 36 in / 122 x 91 cm
Handrail		●		●
Power	110V – 240V AC universal switching power supply (included)			
Shipping Weight	255 lb / 115 kg	265 lb / 120 kg	311 lb / 141 kg	321 lb / 146 kg
UPC Code	809161203707	809161203806	809161203905	809161204001

SOLACE SERIES IN-FLOOR DIALYSIS SCALES



Model	ID-4848S-855RMP	ID-4848SH-855RMP	ID-7248S-855RMP	ID-7248SH-855RMP
Capacity	1,000 lb x .2 lb / 450 kg x .1 kg			
Platform Size	48 x 48 in / 122 x 122 cm	48 x 48 in / 122 x 122 cm	72 x 48 in / 183 x 122 cm	72 x 48 in / 183 x 122 cm
Handrail		●		●
Power	110V – 240V AC universal switching power supply (included)			
Shipping Weight	388 lb / 176 kg	398 lb / 181 kg	512 lb / 232 kg	522 lb / 237 kg
UPC Code	809161204100	809161204209	809161204308	809161204407

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.



203 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153

Mobile: m.detecto.com

SOLD BY: