



# Wheelchair Scale

## Model DS9100

The ultimate in versatility, the heavy-duty DS9100 can serve as a wheelchair scale as well as a handrail or stand-on scale – all in a compact, portable design. The 1,000 lb capacity scale platform is large enough to accommodate any wheelchair, including bariatric wheelchairs. When the patient is on the scale, the ramp does not touch the floor, ensuring that an accurate weight is displayed.

### MAXIMUM CAPACITY

1000 x 0.1 lb (454 x 0.05 kg)

### ACCURACY

+/- 0.1 lb (0.05 kg)

### PLATFORM SIZE

33" W x 32.5" D x 3" H

### SCALE DIMENSIONS - RAMP EXTENDED

49" W x 42" D x 45.5" H

### DISPLAY

1" High LCD

### OPTIONS

RS-232 or USB Communications,  
Second Ramp, and AC Adapter

### BMI CALCULATOR

Input data in FT/IN or CM

### POWER

Six AA Batteries (Included)

### BATTERY LIFE

Two Years


### CONSTRUCTION

Aluminum and Mild Steel Construction  
with Plastic Trim

### WARRANTY

Three Years

The wheelchair tare memory function can store and recall the weights of up to 10 different wheelchairs, eliminating hand calculations to determine the patient's weight. In addition, you can simply enter the tare weight of wheelchairs through the numeric keypad.

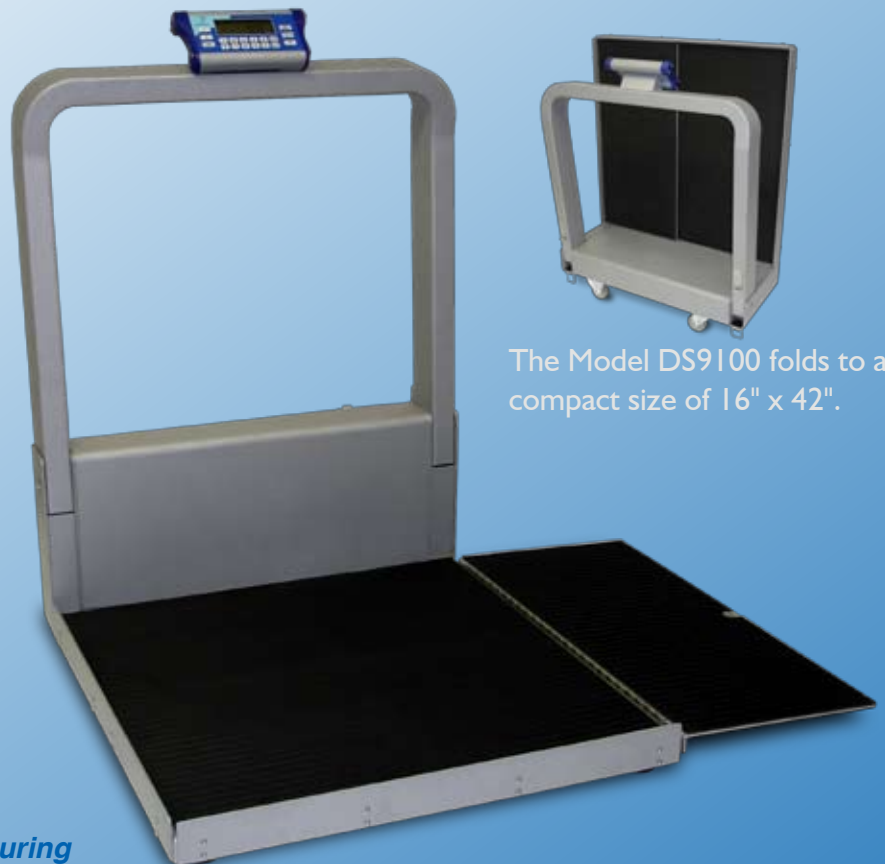

The Model DS9100 folds to a compact size of 16" x 42".



**Featuring**  
**Motion Lock**  
**Technology**

Entering up to 10 wheelchair  
tare weights is easy with the  
numeric keypad

