

NEW & IMPROVED "U" SHAPE FRAME



10200-4

Deluxe Folding Walker, Two Button

10200-4

Adult, 4/cs

10201-2*

Junior, 2/cs

10201-4

Junior, 4/cs

350 LB

WEIGHT CAPACITY

FEATURES:

- Easy push-button mechanisms may be operated by fingers, palms or side of hand.
- Each side operates independently to allow easy movement through narrow spaces and greater stability while standing.
- Unique U-shape frame design creates greater clearance. The patient is able to bring the walker closer to assist in standing.
- Sturdy 1" diameter anodized, extruded aluminum construction ensures maximum strength while remaining lightweight.
- Vinyl contoured hand grip for added comfort.
- Limited Lifetime Warranty.



NEW & IMPROVED FRONT CROSS FRAME!

NEW & IMPROVED UNIVERSAL MODEL!

NEW EASY RELEASE CLOSURE!



300 LB

WEIGHT CAPACITY

UNIVERSAL ADJUSTS FROM 28" - 38"

SIDE BAR ADJUSTS WITH HEIGHT TO PROVIDE ADDITIONAL SUPPORT.



10248N-1

(For 10248N)

Universal (Adult/Junior) Deluxe Folding Walker, Two Button

10248N-1

Universal, 1/cs, For use as Adult and Junior Size, 1/cs, (Adjustable range from 28" - 38")

10248N-2*

Universal, 2/cs, For use as Adult and Junior Size (Adjustable range from 28" - 38")

10248N-4

Universal, 4/cs, For use as Adult and Junior Size, (Adjustable range from 28" - 38")




THIS ITEM IS AVAILABLE IN A 1-PACK!

FEATURES:

- Walker height can be adjusted from 28" - 38". Two walkers in one.
- Side brace adjusts with height of walker to provide additional stability.
- Easy push-button mechanisms may be operated by fingers, palms or side of hand.
- Each side operates independently to allow easy movement through narrow spaces and greater stability while standing.
- Sturdy 1" diameter anodized, extruded aluminum construction ensures maximum strength while remaining lightweight.
- Vinyl contoured hand grip for added comfort.
- Limited Lifetime Warranty.

READ USER'S MANUAL THOROUGHLY BEFORE ADJUSTING.

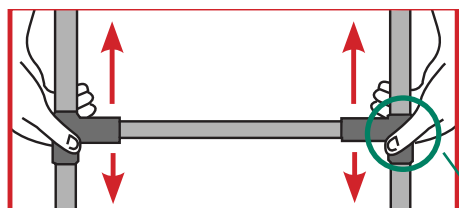

NEW EASY RELEASE CLOSURE!



CAUTION
SLIDE SIDE BRACE USING TWO HANDS ON BOTH ENDS WITH EQUAL FORCE SIMULTANEOUSLY. SIDE BRACE SHOULD REST ON TOP OF PLASTIC SILENCER.

SPECIFICATIONS

Item #	Width (inside hand grip)	Width (inside back legs)	Overall Width	Depth (opened at base)	Height (highest setting)	Height (lowest setting)	Product Weight	Carton Shipping Weight	Weight Limit
10200-4	17"	19"	21.75"	16.5"	39"	32"	6.4 lbs.	28.1 lbs.	350 lbs.
10201-2/4	17"	19"	21.75"	16.5"	32"	25"	6.2 lbs.	14.9/27.8 lbs.	350 lbs.
10248N-1/2	17"	19"	22"	17"	38"	28"	6.5 lbs.	7.5/15 lbs.	300 lbs.
10248N-4	17"	19"	22"	17"	38"	28"	6.5 lbs.	28.5 lbs.	300 lbs.

Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.

walkers

drive

WWW.DRIVEMEDICAL.COM

To order call toll free: 877.224.0946