


Top Layer of Visco Foam with
Dedicated Die-cut Head Section
for Superior Immersion Therapy

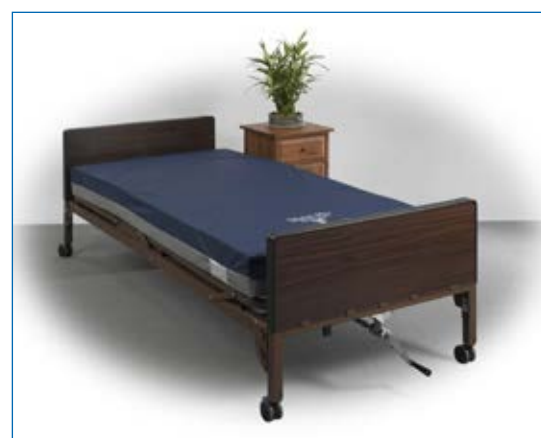

Multi-Ply ShearCare1100 Pressure Redistribution Foam Mattress

Item #	Product Description	Product Dimensions	Product HCPCS	UOM
1100SC-1-FB	Multi-layered/Multi-zoned Foam Mattress	76" (L) x 36" (W) x 6" (H)	E0184	1/ea
1100SC-2-FB	Multi-layered/Multi-zoned Foam Mattress	80" (L) x 36" (W) x 6" (H)	E0184	1/ea
1100SC-3-FB	Multi-layered/Multi-zoned Foam Mattress	84" (L) x 36" (W) x 6" (H)	E0184	1/ea
1100SC-1-RR-FB	Multi-layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out	76" (L) x 36" (W) x 6" (H)	E0184	1/ea
1100SC-2-RR-FB	Multi-layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out	80" (L) x 36" (W) x 6" (H)	E0184	1/ea
1100SC-3-RR-FB	Multi-layered/Multi-zoned Foam Mattress w/3" Elevated Perimeter and Cut-out	84" (L) x 36" (W) x 6" (H)	E0184	1/ea

- Multi-layered pressure redistribution foam mattress with a top layer of Memory-III™ Visco Elastic, an all virgin self-molding and shock-absorbing foam which contours to the body offering reduced interface pressure readings.
- The unique Uni-heel™ section, is a one-piece die-cut high-resilient, low ILD foam offering maximum pressure relief in the vulnerable heel section.
- This combination creates a superior pressure redistributing sleep system.
- Bottom layer of high density foam provides a firm foundation.
- Unique SSPT Rail System™ provides solid-edge support, facilitating safe patient transfers.
- Designated head and foot sections prevent installation errors.
- Top Cover is a 2-way stretch nylon which reduces heat buildup and provides a low-shear surface that is anti-microbial and water-resistant
- Concealed zipper and a Barrier Stop Over-Flap prevents liquids from contaminating mattress core.
- Meets Federal Fire Code CFR 16 part 1633.

specifications

Cover Material	Fluid-Resistant Nylon Top Waterproof Vinyl Bottom
Mattress Layers of Foam	3
Side Rails	Yes
Dedicated Heel Section	Yes
Fire Rating	CFR 16 Part 1633 Federal Fire Code
Weight Capacity	300 lbs.
Warranty	7 Year Limited Foam 1 Year Cover