

DELTA™**Ultra Light 1000, Full-Electric Bed**

**New &
Improved
Game Changer**



**BED, SIDE RAILS AND MATTRESS
SAFETY-TESTED BY SGS,
A NATIONALLY RECOGNIZED TESTING
FACILITY, SATISFY ALL APPLICABLE FDA
ENTRAPMENT GUIDELINES.**

15033

beds


To order call toll free: 877.224.0946

drive™


Truly universal with most
head and footboards by
other manufacturers.
This makes managing
inventory much easier!




Head-
board taller than foot-
board. Better fits home
decor!

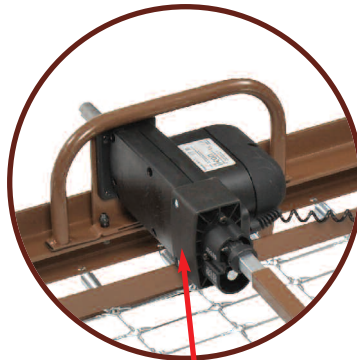

Transition box is on the
foot section NOT on the head
and footboards. That's key to
its interchangeability!

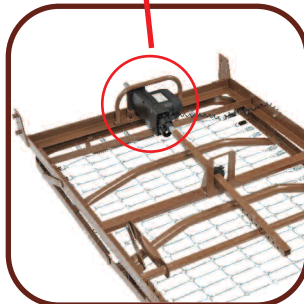

Improved rail placement.
Labels and color-coded
springs assure accurate
installation and patient safety.

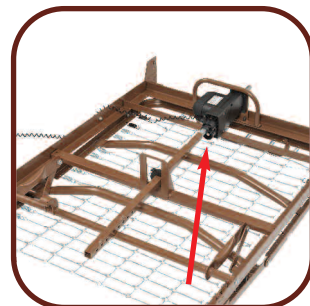

High-low shaft mounted on
bed. Once installed, it can't
be lost or forgotten on fu-
ture installations!

Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.

15033	Ultra Light Full Electric Bed
15033-PKG-11	Ultra Light Full Electric Bed with Full Length Side Rails (item #15001ABV) and 80" Quick 'N Easy Comfort Mattress (item #15076)
15033-PKG-12	Ultra Light Full Electric Bed with Half Rails (item #15201BV) and 80" Quick 'N Easy Comfort Mattress (item #15076)
15033-220	Ultra Light Full Electric Bed, 220 V
15033BV-FR	Ultra Light Full Electric Bed with Full Length Side Rails (item #15001ABV)
15033BV-HR	Ultra Light Full Electric Bed with Half Rails (item #15201BV)
15033BV-PKG	Ultra Light Full Electric Bed with Full Length Side Rails (item #15001ABV) and 80" Inner Spring Mattress (item #15006)
15033BV-PKG-T	Ultra Light Full Electric Bed with Full Length Side Rails (item #15001ABV) and 80" Therapeutic Support Mattress (item #15019)
15033BV-PKG-1	Ultra Light Full Electric Bed with Half Rails (item #15201BV) and 80" Inner Spring Mattress (item #15006)
15033BV-PKG-1-T	Ultra Light Full Electric Bed with Half Rails (item #15201BV) and 80" Therapeutic Support Mattress (item #15019)
15033BV-PKG-2	Ultra Light Full Electric Bed with Full Length Side Rails (item #15001ABV) and 80" Foam Mattress (item #15007)
15096-SEMI	Ultra Light Semi-Electric Bed Reclamation Kit, 1/ea (Includes sleep surface & motor only, no head and foot boards)
15033NKIT	Delta Conversion Kit, 1/ea (Converts Delta Semi Electric bed to Delta Full Electric bed. Includes 2 motors, hand control and hardware.)

FEATURES:

- The Delta 1000 is truly universal because the headboards and footboards are completely interchangeable with Drive's and most other manufacturers bed ends, old and new.
- The transition box mounted on the foot section is key to interchangeability.
- The transition box allows you to change the rotation of the high-low shaft.
- Headboard is taller than footboard to better fit home decor.
- The wood grain panels are break and scratch resistant.
- Foot section weighs 50% less than conventional foot sections to make deliveries safe and easy.
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations.
- Labels and color-coded springs assure accurate installation of side rails for patient safety.
- Unique motor is completely self-contained to reduce weight and noise.
- The motor assembly can be installed and removed with patient in bed.
- The motor is UL approved.
- In the event of a power failure, a 9 volt battery, located in the motor, can lower the head and foot sections nine times. No crank needed.
- New and improved hand pendant has large, easy to use controls.
- Bed ships in two cartons.
- Channel frame construction provides superior strength and reduced weight.
- Zinc coated spring deck.
- Extension kit available model # (15030EXTKIT) extends length from 80" to 84".
- **Warranty: (Original purchaser)**
5 years on motor and hand control, Limited Lifetime on welds and frame,
1 year on all other parts and components.

SPECIFICATIONS

Item #	Bed Ht Low Inside Caster Socket	Bed Ht High Inside Caster Socket	Bed Ht Low Outside Caster Socket	Bed Ht High Outside Caster Socket	Overall Length	Overall Width Weight	Product Weight	Carton Shipping	Weight Limit
15033	9.75"	20"	22.75"	15.25"	88"	36"	162 lbs.	175 lbs.	450 lbs.

Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.