

## GRAVITY 8 DELUXE LONG TERM CARE PRESSURE REDISTRIBUTION MATTRESS

The Gravity 8 multi-layered foam system provides outstanding zoned pressure redistribution, with an emphasis on managing patients at risk for heel pressure ulcers. With a gently sloping heel section, and channel cuts in the middle layer of foam under the foot section, heel pressure off-loading and comfort are maximized.

- 15870 80" (L) x 36" (W) x 6" (H), 1/ea
- 15888 80" (L) x 36" (W) x 6" (H) with Raised Side Rails, 1/ea
- 15884 84" (L) x 36" (W) x 6" (H), 1/ea
- 15885 84" (L) x 36" (W) x 6" (H) with Raised Side Rails, 1/ea
- 15886 76" (L) x 36" (W) x 6" (H), 1/ea
- 15887 76" (L) x 36" (W) x 6" (H) with Raised Side Rails, 1/ea

**CUSTOM SIZES AVAILABLE**

**HCPCS code: E0184**

### COMPARE TO:

\*Direct Supply® Panacea®  
Invacare® Therapy 2000  
Mason® 6500  
Span America® GeoMatt Max  
Blue Chip® Elite  
Joerns® Prevamatt Deluxe

### FEATURES:



- ◆ Outstanding combination of multi-layered die-cut and memory foam available with or without raised side rails.
- ◆ Mattress includes a 6" to 5" heel slope and a "channel cut" center layer of foam under the foot section to target heel ulcer prevention.
- ◆ Top layer uses unique die-cutting and a true memory foam heel section to optimize zoned pressure redistribution.
- ◆ Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress significantly.
- ◆ Cover is fire retardant, fluid proof, tear resistant, low sheer, and easy to clean.
- ◆ Latex free.
- ◆ 350 lb. Weight Capacity.
- ◆ **1 Year Warranty on the cover.**
- ◆ **10 Year Limited Warranty.**

**\*\*Always use mattress sizes that correspond to the exact mattress deck surface size, and secure with proper mattress stops.**

\*Direct Supply® and Panacea® are registered trademarks of Direct Supply, Inc.


**NEW!**

**LATEX FREE!**

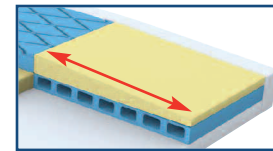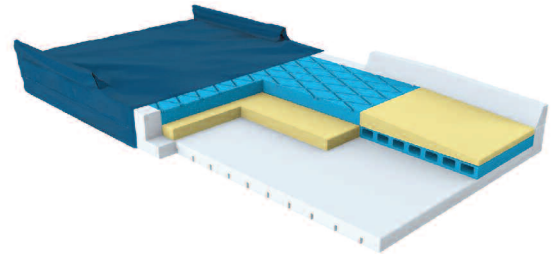

**HEEL SLOPE  
AND CHANNEL CROSS CUTS**

**PREVENT HEEL PRESSURE ULCERS  
WITH THE GRAVITY 8'S HEEL SLOPE,  
CHANNEL CUTS, AND OPEN-CELL  
VISCO-ELASTIC MEMORY FOAM  
FOOT SECTION!**

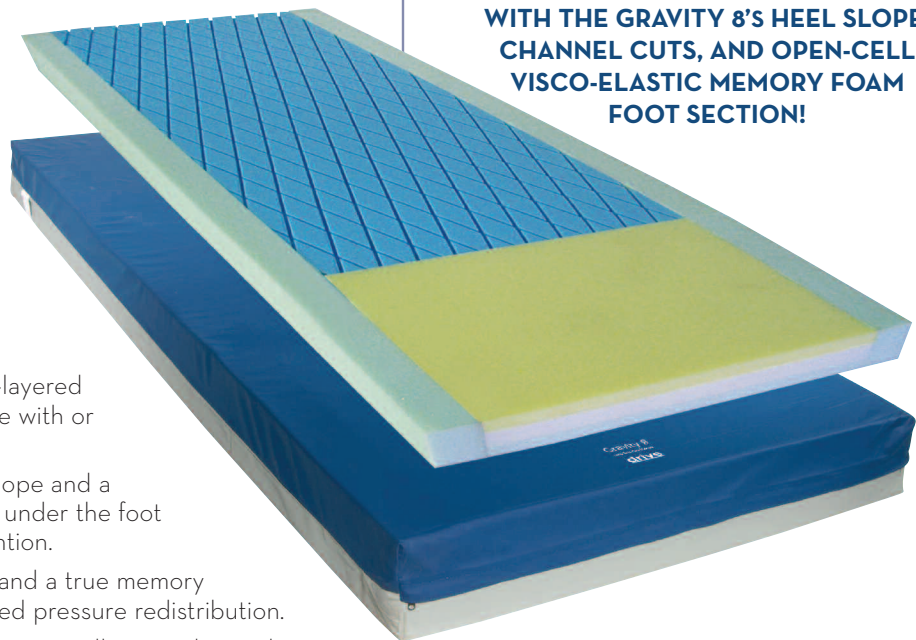

**To order, call LTCCare Customer Service Toll Free at (877) 374-8763,  
contact your Sales Representative or visit us at [www.drivemedical.com](http://www.drivemedical.com) TODAY!**