


# IMAGE EC



Shown with  
ECEL-R-P  
Seat Plate.  
Sold Separately.

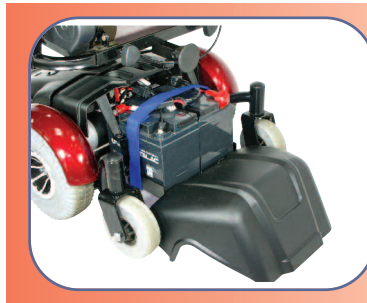

EASY TO  
ACCESS  
BATTERIES!



## Image EC Mid Wheel Drive


Available  
**1 Pack**

2800ECBL-RCL  
2800ECBU-RCL  
2800ECBL-RCL-20  
2800ECBU-RCL-20

18", Blue, (Replaces 2800GTBL-RCL), 1/cs  
18", Burgundy, (Replaces 2800GTBU-RCL), 1/cs  
20", Blue, 1/cs  
20", Burgundy, 1/cs

### FEATURES:

- Mid wheel drive allows for tighter turning radius.
- PG Controller with diagnostic fault codes.
- Off board charger.
- Joystick can be mounted on either side.
- Front swivel anti-tippers have true independent adjustable suspension. Helps prevent the drive wheels from losing contact with the ground.
- Drive mid wheels are flat free including coil spring tension suspension.
- Uniquely designed Dynamic Drive suspension for greater climbing ability and comfort when traveling forward and in reverse.
- Front and Rear semi-pneumatic swivel wheels.
- Foot plate is angle, height and depth adjustable.
- Elevating legrest or footrest kit available (SFELR-KIT).
- 18" deluxe height adjustable swivel captain's seat. (20" optional)
- Positioning belt standard.
- Arms are width, height and angle adjustable.
- Standard headrest is height adjustable.
- Seat reclines for comfort and folds down for portability.
- Has been tested and surpassed the ANSI RESNA testing requirement.
- Comes with batteries (Item #U1)
- Batteries are accessible from the rear of the chair.
- 300 lb. Weight Capacity.

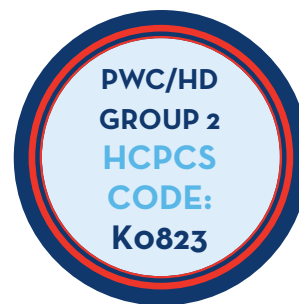

Blue



Burgundy

\*Important Notice: Item numbers that display an asterisk (\*) reflect an item that will be discontinued or package change.