

## Technical data

TECHNICAL DATA AquaJoy Premier Plus & AquaJoy Junior Plus  
Width of frame 385mm (15.2 inches) Length of seat: 580 mm (22.8 inches) Width of seat with extended seat side flaps: 686 mm (27 inches) Width of seat with seat side flaps folded in: 390 mm (15.3 inches) Height of seat in its lowest position: 80 mm (3.1 inches) Height of seat in its highest position: 460 mm (18.1 inches) Depth of seat measured from Backrest to front seat edge: 490 mm (19.3 inches) Height of backrest: 604 mm (23.8 inches) Width of backrest: 340 mm (13.4 inches) Carrying capacity: 26.7 stone (170 kgs) Total Weight without battery: 12.5kg (27.5 lb) Weight of seat: 7.4kg (16.3 lb) Weight of backrest: 5.1kg (11.2 lb) Weight of battery handcontroller: 664g Duration of stroke (80 kg load): Up approx. 50sec, Down approx. 35sec Operating voltage: 13.2V DC Maximum current consumption: 6A Charging rate: Nominal Input: 13.2V DC Motor noise with no load: approx. 50dB (according to DIN 3746) Medical equipment type B Class Protection: IP 66 Duty cycle: Max. 5%, Max. 3min./hour Battery: 13.2V 1300mAh (Ni-MH) Battery Charger Input: 240V AC/50Hz, 230V AC/60Hz, 110V AC/60Hz Output: 18V DC 250mA/4.5VA Weight: 0.3kg(UK 282g /EU 285g /USA 275g) Charging time with flat battery: approx. 12h



**Warranty:** The AquaJoy Premier Plus and AquaJoy Junior Plus come with a three year guarantee excluding battery hand control and covermat (1 year). (Unless otherwise negotiated)

## Premier Plus & Junior Plus Accessories



### ✓ Swivel & Slide:

Used to assist the transfer on and off the bathlift.



### ✓ Pommel:

Prevents slipping.



### ✓ Wedge Cushion:

Sits under knees to prevent slipping.



### ✓ Chest & Lap Harness:

Ensuring safety at all times.



### ✓ Headrest:

Provides extra support and comfort whilst bathing.



### ✓ Upper Body Wing Support:

Provides extra support and security for the upper body.



### ✓ Unique Gel Filled Covermats:

Provides extra padding and comfort.



### ✓ Carrier Bag:

The ideal and easy way to store your bathlift. Or if you go on holiday, why not take your bathlift with you!



### ✓ Middle Suction Cups:

Provides increased suction for baths with dimples.



AquaJoy Bathlifts Ltd  
Consett Business Park  
Consett  
County Durham  
DH8 6BN  
United Kingdom

Tel: +44 (0)1207 501555

Fax: +44 (0)1207 599789

email: [info@aqua-joy.com](mailto:info@aqua-joy.com)

[www.aqua-joy.com](http://www.aqua-joy.com)



# Why choose any other bathlift?

"Amazing, particularly cost!"

"At last, a bathlift that does everything!"

"A really good product which is very user friendly."

"I would not hesitate to recommend this versatile product."

"Just what we have been waiting for, an all singing, all dancing bathlift"

"AquaJoy has designed a product that benefits everyone involved"

"Clients find the hand control extremely easy to use and understand"

"An excellent and versatile product."

"I can't believe how much legroom clients have with this bathlift – its fantastic"

"Your bathlifts are ultra reliable. Thank you for a refreshingly great service."

"An excellent product - ideally priced for social services and private clients"

"In this day and age, to receive such high levels of service is rare - thankyou AquaJoy"

"I am really impressed with everything! AquaJoy have listened!"

"The best bathlift I have seen to date"

**NOW WITH NEW ENHANCED FEATURES**

AquaJoy Bathlifts Ltd – AN AWARD WINNING COMPANY

AquaJoy Premier Plus & Junior Plus Bathlifts



# THE AQUAJOY PREMIER PLUS AND JUNIOR PLUS BATHLIFTS – Developed for you, by you

AquaJoy Bathlifts Ltd has established itself as a brand you can rely on within the Bathlift market place. The award winning company provides high quality, versatile, reliable products supported by unbeatable levels of customer service. The AquaJoy Premier Plus & Junior Plus Bathlifts were developed and designed after close consultation with healthcare professionals, loan store managers and end users. Because of this, the Bathlifts have become firm favourites throughout the UK and Europe.

The AquaJoy Premier Plus & Junior Plus Bathlifts fit into almost any style of bath, maximise legroom and are easy to put in and take out of your bath without requiring any adaptation or modification to your bathroom.

AQUAJOY JUNIOR PLUS

## So why choose any other bathlift?

### SUITABLE FOR ADULTS AND CHILDREN

- Adaptable accessories for users with varying ages and needs

### BARIATRIC WEIGHT LIMIT AS STANDARD 170KG (26.7 STONE)

- Ideal for heavier clients

### RECLINING BACKREST AT THE TOUCH OF A BUTTON



- Allows the client to recline in the bath like they used to

### FITS ALMOST ANY STYLE OF BATH

- Suitable for deep baths – maximum seat height 460mm
- Suitable for shallow baths – minimum seat height 80mm

### MAXIMISES LEGROOM IN THE BATH

- Bathlift seat locates as far back in the bath as possible
- Recline angle makes full use of all bath space

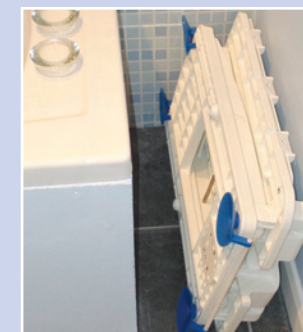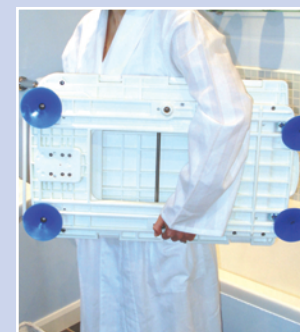

## Moving and Handling – made easy!

### TWO PIECE RECLINING BATHLIFT

- Separates easily into two compact sections
- No awkward locks or buttons to press, backrest simply folds forward to detach
- Compact design saves storage space and allows the Bathlift to be carried against the body for an unobstructed view

### EASY RELEASE SUCTION FEET WITH LIFTING BAR AS STANDARD – NOW BLUE!

- Limits stooping and stretching when lifting the Bathlift in and out of the bath
- No struggling under the bathlift seat to detach individual suction feet

## Easy to clean and maintain

### CAN BE USED WITH OR WITHOUT COVERMATS

- Detachable and can be washed in a standard washing machine up to 60 degrees protecting against unwanted germs and bacteria

### SMOOTH, CREVICE FREE SURFACE

- Dirt and used water do not get trapped
- Makes cleaning quick and easy

### STREAMLINED FRAMEWORK

- Provides access to the Bathlift and bath for cleaning when in situ



## Simple and safe to operate

### LIGHTWEIGHT, FULLY WATERPROOF, FLOATING HANDCONTROL

- Comfortable to hold and operate
- Suction pads on reverse to stick to the side of the bath or to tiles for ease of use

### CONVEX AND CONCAVE BUTTONS

- Soft to touch
- Easy to operate for users with limited dexterity or visual impairment

### SIMPLE CHARGING PROCESS

- New and improved pronged connectors for ease of use

### BATTERY INDICATOR LIGHT

- Indicates when recharging is required

### NEW CHARGER INDICATOR LIGHT

- Indicates charging is in progress

### SAFETY FEATURE

- The Bathlift will not lower unless there is sufficient charge to raise the client back to the top of the bath

AQUAJOY PREMIER PLUS

“At last, a Bathlift that does everything!” “AquaJoy has designed a product that benefits everyone involved” “An excellent and versatile product” “The best Bathlift I have seen to date!”