

# Adult Anterior Safety Roller

- CE 1000 B Adult, Burgundy, 1/cs
- CE 1000 BK Adult, Black, 1/cs
- CE 1000 NBL Adult, Narrow Width, Blue, 1/cs

## FEATURES:

- Applying pressure to the handlebar activates braking mechanism. When brake is engaged, the brake stopper hits the wheel and the wheels cannot roll. Ideal for patients with limited hand function and/or limited cognizance. The Safety Roller halts in step with the patient and the braking mechanism automatically engages if patient stumbles.
- Height adjustable in 1" increments.
- Adjustable brake spring tension.
- Patient is situated in the center of the Safety Roller, which provides maximum support and improved posture.
- Skid resistant, non-marring, easy rolling neoprene wheels. Ideal for indoor and outdoor use.
- Folds easily and stands on its own in the folded position.
- Limited Lifetime Warranty.



CE 1000 BK

**400 LB. WEIGHT CAPACITY!**


**FOLDS EASILY!**



ITEM #	HANDLEBAR HEIGHT RANGE	OVERALL WIDTH	WIDTH INSIDE HANDLE GRIP	DEPTH OPENED AT BASE	PRODUCT WEIGHT	CARTON SHIPPING WEIGHT	WEIGHT LIMIT
CE 1000 B/BK	29" - 36"	21.5"	17"	25"	16 lbs.	21 lbs.	400 lbs.
CE 1000 NBL	28" - 35"	20"	15"	25"	12 lbs.	17 lbs.	400 lbs.

\* For a complete line of accessories, see pages 426-430.

\*Important Notice: Item numbers that display an asterisk (\*) reflect an item that will be discontinued or package change.

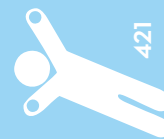

421

mobility

safety rollers

adult

**Wenzel**iteRe/hab

To order call toll free: 877.224.0946

**drive**

WWW.DRIVEMEDICAL.COM