

# Otter Bathing System



**Available**  
**1 Pack** Model #: Sold 1 per case



## Standard Fabric

OT 1000  
OT 2000  
OT 2000

Small  
Medium  
Large

## Standard Fabric

OT 1000 SF  
OT 2000 SF  
OT 2000 SF

Small  
Medium  
Large

**Weight Capacity:**

60/120/160 lbs

**Warranty:**

Limited Lifetime


Accessories see page 48

## Features

- Plastic frame.
- Fabric can be removed and machine-washed.
- Seat fabric can be slipped off to bathe a child sitting on the floor of the tub, supported by the bath chair back.
- Seat and back are angle adjustable and have five angle adjustments at 0°, 22.2°, 45°, 67.5° & 90°.
- Adjustable slip-resistant legs raise chair up to 7". Height can be adjusted with child in chair.
- Uni-Bars allow seat and back angle adjustments to be made with one hand.
- Leg straps control abduction and adduction.
- Positioning straps can be placed anywhere along the length of the frame and can be removed if not required.
- Head cushions are height and width adjustable and provide support and comfort. They can also be used as lateral supports.
- Adjustments can be made with user in bath chair.
- Folds flat for easy storage.
- Two Year Limited Warranty.



UNI-BAR

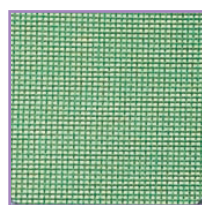

**Standard Fabric**  
Vinyl covered nylon  
dries easily!



**Soft Fabric**  
Polyester knit with a trunk  
strap that assists in lateral  
support. It has extra padding  
between the fabric & the  
Otter frame.

## Specifications

Item #	Seat Depth	Seat Width	Back Height	Seat Height	Overall Outside Width	Overall Outside Length	User Height	Product Weight	Carton Shipping Weight	Weight Limit
OT 1000/SF	13"	14.5"	18"	2" - 7"	17"	36.5"	Up to 36"	11 lbs.	14 lbs.	60 lbs.
OT 2000/SF	13"	14.5"	25"	2" - 7"	17"	42"	32" - 50"	12 lbs.	15 lbs.	120 lbs.
OT 3000/SF	16.5"	14.5"	32"	2" - 7"	17"	50"	46" - 68"	14 lbs.	18 lbs.	160 lbs.

\*Important Notice: Item numbers that display an asterisk (\*) reflect an item that will be discontinued or package change.

To order call toll free: 877.224.0946

WENZELITE


Wenzel's Re/Hub  
**drive**

# Otter Accessories




To order call toll free: 877.224.0946

WENZELITE

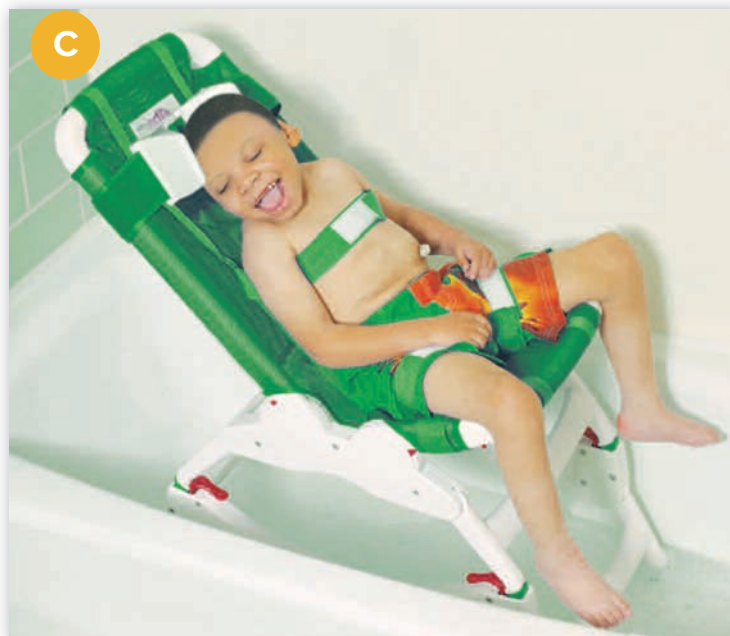

One stand for all three Otter sizes!  
The tub stand raises the Otter height to a convenient height for easier transfers.

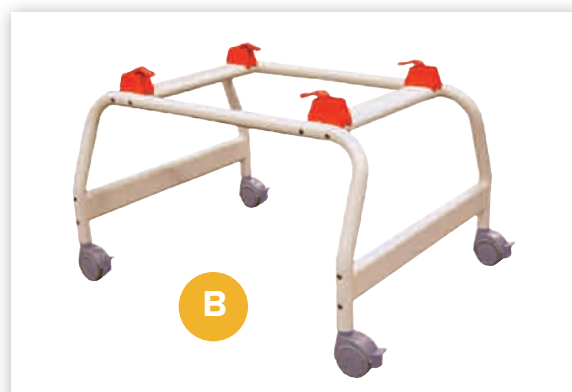

## A Tub Stand OT 8010 1/cs

- Slip-resistant legs.
- Legs fold flat for storage.
- Quick-release locks allow chair to be easily attached or removed.
- Height adjustable (5"-10").
- 200 lb. weight capacity.

## B Shower Stand OT 8020 1/cs

- Sturdy aluminum stand.
- Swivel locking casters.
- Quick-release locks allow chair to be easily attached or removed.
- Ships unassembled.

## C Lateral Supports OT 8000 1/cs

- Wedges and adjustable strap provide lateral support.

## Specifications

Item #	Width	Height	Length	Product Weight	Carton Shipping Weight	Weight Limit
OT 8010	17"	5"-10"	32"	8 lbs.	11 lbs.	200 lbs.
OT 8020	17"	21"	33"	17 lbs.	20 lbs.	200 lbs.

\*Important Notice: Item numbers that display an asterisk (\*) reflect an item that will be discontinued or package change.