

Gait trainers are designed for individuals who do not have the functional gait skills or postural control to safely and independently use a standard walker.

The Benefits of Gait Trainers

- Encourages early cognitive and learning development.
- Cardiopulmonary strengthening and conditioning.
- Reduces flexion contractures.
- Increases spatial awareness and head control to achieve independent ambulation.
- Teaches individual to take steps.

The Wenzelite gait trainers offer varied levels of support and each accessory can be removed or added according to the individual user's needs and abilities. The posterior and anterior models are each available in three sizes.



Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.




introduction

Wenzelite Re/hab

To order call toll free: 877.224.0946

drive

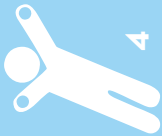

NEW!

Trekker

Trekker Gait Trainer

TK 1000 Tyke, Orange, 1/bx
TK 2000 Junior, Green, 1/bx
TK 3000 Youth, Blue, 1/bx

FEATURES:

- Height adjustable in 1" increments.
- Can be used in the anterior and posterior position.
- Swivel casters that can be locked to non-swivel.
- Variable resistance to control the speed of the wheel's rotation.
- Brakes to lock the wheels to prevent the unit from moving.
- Directional lever allows both forward and reverse mobility or only forward mobility.
- Accessories mount tool-free onto the frame and can be removed or added according to the individual user's needs and abilities.
- Junior & Youth Trekkers fold easily.
- Limited Lifetime Warranty



TK 1000



TK 2000



TK 3000

ITEM #	ADJUSTABLE HEIGHT	OVERALL WIDTH	WIDTH BETWEEN HANDLEBARS	OVERALL LENGTH	CARTON SHIPPING WEIGHT	WEIGHT LIMIT
TK 1000	15" - 18.5"	22.25"	12.5"	32"	29"	75 lbs.
TK 2000	21" - 26"	24.5"	17"	31"	34"	150 lbs.
TK 3000	29" - 39"	26.5"	18"	37"	36"	200 lbs.

mobility

gait trainers

trekker

Wenzel[®]ite Re/hab

To order call toll free: 877.224.0946

drive

ACCESSORIES:

Forearm Platform

TK 1035 S Forearm Platforms, Small, 1/bx

TK 1035 L Forearm Platforms, Large, 1/bx

FEATURES:

- Height & depth adjustable.
- Platforms can rotate, supinate and pronate.
- Handgrips extend forward and backward, supinate, pronate and are angle adjustable.
- Can be mounted anywhere on the handlebar.
- Can be mounted facing toward the inside or outside of the handlebar for lateral adjustability.
- Contoured armrest with arm and wrist straps.

Handgrips

TK 1045 Handgrips, 1pr/bx

Handgrips can be mounted anywhere on the handlebar and they are height adjustable and rotate.

Ankle Prompts

TK 1060 S Ankle Prompts for TK 1000, 1pr/bx


TK 1060 M Ankle Prompts for TK 2000, 1pr/bx

TK 1060 L Ankle Prompts for TK 3000, 1pr/bx

FEATURES:

- Used for patients with spasticity, poor motor control, ataxia and excessive adduction.
- Length adjustable, padded ankle straps maintain legs and/or feet in abduction and prevent scissoring.
- Spring clips control stride length.

NEW!

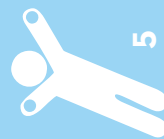

mobility

gait trainers

trekker

Wenzel *ite Re/hab*

To order call toll free: 877.224.0946

drive



Hip Positioner & Pad

TK 1070 S Hip Positioner & Pad, Small, 1/bx
TK 1070 L Hip Positioner & Pad, Medium, 1/bx



FEATURES:

- Height adjustable seat harness positions hip and pelvis and allows user to shift weight from one leg to the other.
- Shorten length of straps to assist in weight bearing or lengthen straps to use as a safety sling.
- Position of straps on handlebar will affect user's pelvis position.
- Removable pad provides cushioned comfort.



Trunk Support

TK 1080 S Trunk Support, Small, 1/bx
TK 1080 M Trunk Support, Medium, 1/bx
TK 1080 L Trunk Support, Large, 1/bx

FEATURES:

- Secures the torso at desired angle.
- Height, depth, width and angle adjustable.
- Padded for cushioned comfort.
- Stays open for easy transfers.
- Can be mounted towards the inside or outside of handlebar.
- Straps adjust the fit and the forward leaning and angle of the user.
- In its lowest position, can be used as pelvic positioner.
- Adjustments can be made while user is in the gait trainer.



Thigh Prompts

TK 1090 S Thigh Prompts, Small, 1pr/bx
TK 1090 L Thigh Prompts, Large, 1pr/bx

FEATURES:

- Height adjustable
- Can be positioned further from or closer to the user to limit adduction/abduction.
- Moves with user's gait pattern allowing weight shifting from one leg to the other.
- Prevents leg scissoring.