

## 3M™ Tegaderm™ Transparent Film Roll

*All the benefits of Tegaderm™ Transparent Film Dressings in a non-sterile roll*

# Tegaderm™ Transparent Film Roll



Simple.  
Dependable.  
Trusted.

### **A non-sterile version of Tegaderm™ Transparent Film Dressings**

Tegaderm™ Roll provides maximum versatility and is the affordable choice for skin protection and waterproof cover dressings.

- The roll version allows for maximum versatility.
- Tegaderm™ Film provides a protective cover over at-risk, including newly healed, intact skin helping to prevent skin breakdown caused by friction and body fluids. Tegaderm™ Film is breathable, letting oxygen in and moisture vapor out, allowing the skin to function normally.
- The perfect solution whenever a waterproof cover dressing is the appropriate choice.



# 3M™ Tegaderm™ Transparent Film Roll

## Unique Features



- Thin film is easily cut to desired shape or size. If needed, corners can be shaped and rounded.
- Conformable – Fits body contours, permits free movement without stress on skin.
- Waterproof, bacterial *and viral*\* barrier, helps to prevent external contamination of the site, as well as helps to protect the patient and healthcare workers.
- Easy stretch and release removal.



Easy application system that can be applied with gloved hands.



Waterproof fixation of primary dressings, such as sterile gauze, foam, alginate, hydrogel, ABD pad, etc.



Fixation of drainage and infusion tubing.



Protection of at risk skin from friction and body fluids.

\*In vitro testing shows that the transparent film of Tegaderm™ and Tegaderm™ HP dressings provides a viral barrier from viruses 27 nm in diameter or larger while the dressing remains intact without leakage.

# Application and Removal

## Application Instructions



**Step 1**

Cut to desired length and round corners or cut to shape for given application.



**Step 2**

Remove the backliner. Hold dressing by side tabs. Center over dressing or area requiring fixation or protection and apply without stretching.



**Step 3**

Gently smooth down the dressing, helping to seal the film. Remove the top liner.

## Removal Instructions



**Step 1**

Grasp one edge of the dressing and gently pull straight out to stretch it and release adhesion, carefully holding onto skin with other hand.



**Step 2**

If necessary, press a small piece of tape onto a corner to start an edge of the dressing. By stretching the dressing, the material expands and releases until the dressing is removed.

## Tegaderm™ Roll can also be used for:

- Protection of intact skin during negative pressure wound therapy (i.e. periwound skin)
- Shower cover over primary dressings such as tape and gauze (may require non adherent layer between primary dressing and Tegaderm™ Roll)
- Securement of transparent dressing edges when lift occurs prior to protocol change date
- Protection of temporary skin markings used for radiation therapy

## Affordability

A lower cost option for many applications. Non-sterile 3M™ Tegaderm™ Roll is an affordable product choice compared to sterile separate dressings, where appropriate. It can be significantly less expensive than traditional transparent film dressings and should be considered as a first choice for skin protection and waterproof cover dressings.




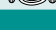









# 3M™ Tegaderm™ Transparent Film Dressing Family

## Ordering Information

CHG

Antimicrobial

NEW!  
Smaller  
Size

Catalog No. NHRIC No.	Size	Dressings/ Box	Boxes/ Case	HCPCS Code*
3M™ Tegaderm™ CHG Chlorhexidine Gluconate IV Securement Dressing				
 <b>1657</b> NHRIC 8333-1657-01	3 1/2 in x 4 1/2 in 8,5 cm x 11,5 cm	25	4	A6257
 <b>1658</b> NHRIC 8333-1658-01	4 in x 4 3/4 in 10 cm x 12 cm	25	4	A6257
 <b>1659</b> NHRIC 8333-1659-01	4 in x 6 1/4 in 10 cm x 15,5 cm	25	4	A6258
 <b>1660</b> NHRIC 8333-1660-01	2 3/4 in x 3 3/4 in 7 cm x 8,5 cm	25	4	A6257
Catalog No. NDC No.	Size	Dressings/ Box	Boxes/ Case	HCPCS Code*
3M™ Tegaderm™ IV Transparent Film Dressing with Border				
 <b>1610</b> NDC 8333-1610-01	2 in x 2 3/4 in 5 cm x 5,7 cm	100	4	A6257
 <b>1650</b> NDC 8333-1650-01	4 in x 6 1/4 in 10 cm x 15,5 cm	25	4	A6258
 <b>1655</b> NDC 8333-1655-01	3 1/2 in x 4 1/2 in 8,5 cm x 11,5 cm	50	4	A6257
 <b>1633</b> NDC 8333-1633-01	2 3/4 in x 3 3/4 in 7 cm x 8,5 cm	100	4	A6257
 <b>1635</b> NDC 8333-1635-01	3 1/2 in x 4 1/4 in 8,5 cm x 10,5 cm	50	4	A6257
3M™ Tegaderm™ Transparent Film Dressing with Border				
 <b>1614</b> NDC 8333-1614-05	2 1/2 in x 2 3/4 in 6 cm x 7 cm	100	4	A6257
 <b>1616</b> NDC 8333-1616-05	4 in x 4 3/4 in 10 cm x 12 cm	50	4	A6258
3M™ Tegaderm™ HP (Holding Power) Transparent Film Dressing Frame Style				
 <b>9519HP</b> NDC 8333-9519-01	2 3/4 in x 2 3/4 in 6 cm x 6 cm	100	4	A6257
Catalog No. NHRIC No.	Size	Dressings/ Box	Boxes/ Case	HCPCS Code*
3M™ Tegaderm™ IV Securement Dressing Designed for the BD Nexiva™ Closed IV Catheter System				
 <b>9525HP</b> NHRIC 8333-9525-01	2 1/2 in x 2 3/4 in 6,5 cm x 7 cm	100	4	A6257

**3M™ Tegaderm™** SIMPLE. DEPENDABLE. TRUSTED.  
Dressings

For more information visit our website at [www.3M.com/tegaderm](http://www.3M.com/tegaderm), contact your 3M Skin & Wound Care representative, or call the 3M Health Care Customer Helpline at **1-800-228-3957**.

These products can be ordered from your local distributor. Outside the United States contact the local 3M subsidiary.

\*HCPCS codes as referenced in the Local Coverage Determination Surgical Dressing Policy.

### Disclaimer:

HCPCS codes have been provided to assist you in the preparation of Medicare Part B claims. Please note, however, that the reimbursement information provided by 3M Health Care and its representatives is intended to provide general information relevant to coverage and coding for 3M products. Insurers' reimbursement policies can vary and the use of the codes discussed here does not guarantee that an insurer will cover or pay at any particular level. Health care providers should exercise independent clinical judgment in choosing the codes which most accurately describe the products provided.























**3M**

**Skin & Wound Care Division**  
**3M Health Care**  
3M Center, Building 275-4W-02  
St. Paul, MN 55144-1000  
USA  
1-800-228-3957  
[www.3M.com/healthcare](http://www.3M.com/healthcare)

**3M Canada Company**  
P.O. Box 5757  
London, Ontario N6A 4T1  
Canada  
1-800-364-3577  
[www.3M.com/ca/healthcare](http://www.3M.com/ca/healthcare)

Nexiva is a trademark of Becton, Dickinson and Company.

3M and Tegaderm are trademarks of 3M. Please recycle. Printed in U.S.A.  
© 3M 2005, 2006, 2007, 2009. All rights reserved.  
70-2009-7006-2 (119.4) DPI

Catalog No. NDC No.	Size	Rolls/ Box	Boxes/ Case	HCPCS Code*
3M™ Tegaderm™ Transparent Film Roll (Non-Sterile)				
 <b>16002</b> NDC 8333-1600-20	2 in x 11 yd 5 cm x 10 m	1	4	A6257-A6259
 <b>16004</b> NDC 8333-1600-40	4 in x 11 yd 10 cm x 10 m	1	4	A6257-A6259
 <b>16006</b> NDC 8333-1600-60	6 in x 11 yd 15 cm x 10 m	1	4	A6257-A6259
Catalog No. NDC No.	Size	Dressings/ Box	Boxes/ Case	HCPCS Code*
3M™ Tegaderm™ Transparent Film Dressing Frame Style				
 <b>1622W</b> NDC 8333-1622-05	1 3/4 in x 1 3/4 in 4,4 cm x 4,4 cm	100	4	A6257
 <b>1624W</b> NDC 8333-1624-05	2 3/4 in x 2 3/4 in 6 cm x 7 cm	100	4	A6257
 <b>1626W</b> NDC 8333-1626-05	4 in x 4 3/4 in 10 cm x 12 cm	50	4	A6258
 <b>1630</b> NDC 8333-1630-05	4 in x 4 1/2 in 10 cm x 11,5 cm	50	4	A6258
3M™ Tegaderm™ Transparent Film Dressing Frame Style (Small quantity packages)				
 <b>9505W</b> NDC 8333-9505-05	2 3/4 in x 2 3/4 in 6 cm x 7 cm	20	10	A6257
 <b>9506W</b> NDC 8333-9506-05	4 in x 4 3/4 in 10 cm x 12 cm	10	10	A6258
3M™ Tegaderm™ HP (Holding Power) Transparent Film Dressing Frame Style				
 <b>9534HP</b> NDC 8333-9534-01	2 3/4 in x 2 3/4 in 6 cm x 7 cm	100	4	A6257
 <b>9536HP</b> NDC 8333-9536-01	4 in x 4 3/4 in 10 cm x 12 cm	50	4	A6258
 <b>9543HP</b> NDC 8333-9543-01	4 1/2 in x 4 3/4 in 11,5 cm x 12 cm	12	4	A6258
 <b>9545HP</b> NDC 8333-9545-01	2 1/2 in x 2 1/2 in 5,4 cm x 6,4 cm	50	6	A6257
 <b>9546HP</b> NDC 8333-9546-01	4 in x 4 1/2 in 10 cm x 11,5 cm	50	4	A6258
 <b>9548HP</b> NDC 8333-9548-01	5 1/2 in x 6 1/2 in 14 cm x 16,5 cm	10	8	A6258
3M™ Tegaderm™ Transparent Film Dressing Frame Style				
 <b>1634</b> NDC 8333-1634-01	2 3/4 in x 2 3/4 in 6 cm x 7 cm	100	4	A6257
 <b>1626</b> NDC 8333-1626-01	4 in x 4 3/4 in 10 cm x 12 cm	50	4	A6258
 <b>1627</b> NDC 8333-1627-01	4 in x 10 in 10 cm x 25 cm	20	4	A6258
 <b>1628</b> NDC 8333-1628-01	6 in x 8 in 15 cm x 20 cm	10	8	A6258
 <b>1629</b> NDC 8333-1629-01	8 in x 12 in 20 cm x 30 cm	10	8	A6259
3M™ Tegaderm™ Transparent Film Dressing First Aid Style				
 <b>1620</b> NDC 8333-1620-01	2 3/4 in x 2 3/4 in 6 cm x 7 cm	100	6	A6257
 <b>1621</b> NDC 8333-1621-01	4 in x 4 3/4 in 10 cm x 12 cm	50	4	A6258

Basic

Standard

Roll