


Foam Tipped Applicator

Acetal handle

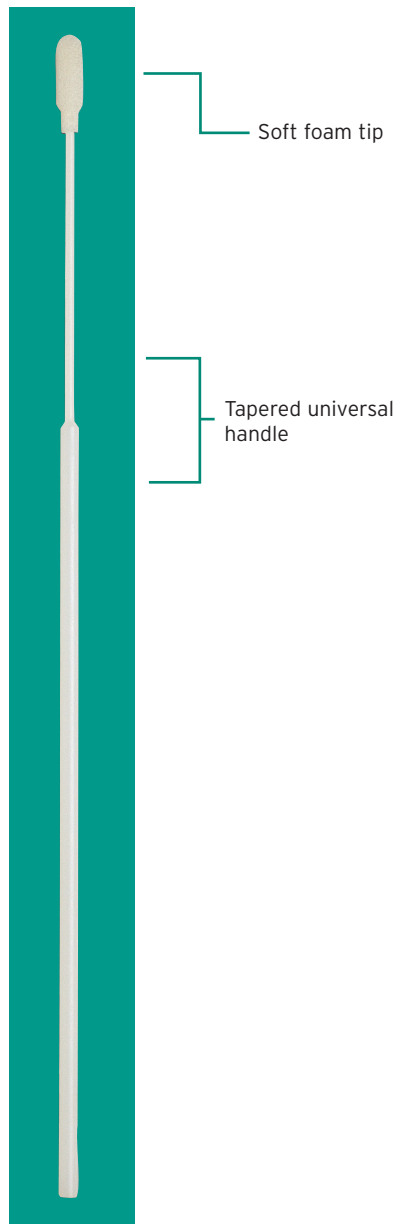

Actual Size

Description

This specialty applicator is tipped with 100ppi medical grade polyurethane foam. It is bonded securely to a tapered acetal handle using a thermoplastic adhesive and heat stake process. Often used as a general swab for medical or veterinary diagnostics.

Dimensions

6.062" Overall Length (153.975 mm)

Tip

- 0.547" Length (13.894 mm)
- 0.141" Diameter (3.581 mm)

Handle

- 5.996" Length (152.298 mm)
tapered to 0.049" (1.244 mm)
- 0.098" Dia. (2.489 mm)

Standards

Non-sterile

- ISO 9001:2008
- FDA GMP's
- ISO 13485:2003
- CE marked for medical use (sterile only)
- Certificate of sterility and/or compliance available with each shipment upon request

Qualities

- Medical grade quality
- Foam has excellent collection and release properties
- Can be used with IPA, DI water, or dry
- Soft small tip for patient comfort
- Lint and residue free

Applications

Industrial / Medical / Veterinary

- Cleaning fiber optic connectors
- Applying or removing adhesives and residues
- Specimen collection
- Diagnostic sampling
- Recovery of micro-organisms and fluids

Ref #	Packaging
1406-PF <i>Non-sterile</i>	20/50 (20 bags per case, 50 swabs per bag)
25-1406 1PF 50 <i>Sterile</i>	10/50/1 (paper/poly pouch) (10 boxes per case, 50 individually wrapped packages per box, 1 swab per package)



Puritan Medical Products Company LLC
Guilford, Maine 04443-0149 USA

Tel 800-321-2313 (US and Canada Only)
207-876-3311

Fax 800-323-4153
207-876-3130

Web www.puritanmedproducts.com

Email sales@puritanmedproducts.com

1-800-321-2313
puritanmedproducts.com