


Foam Tipped Applicator

Polystyrene handle

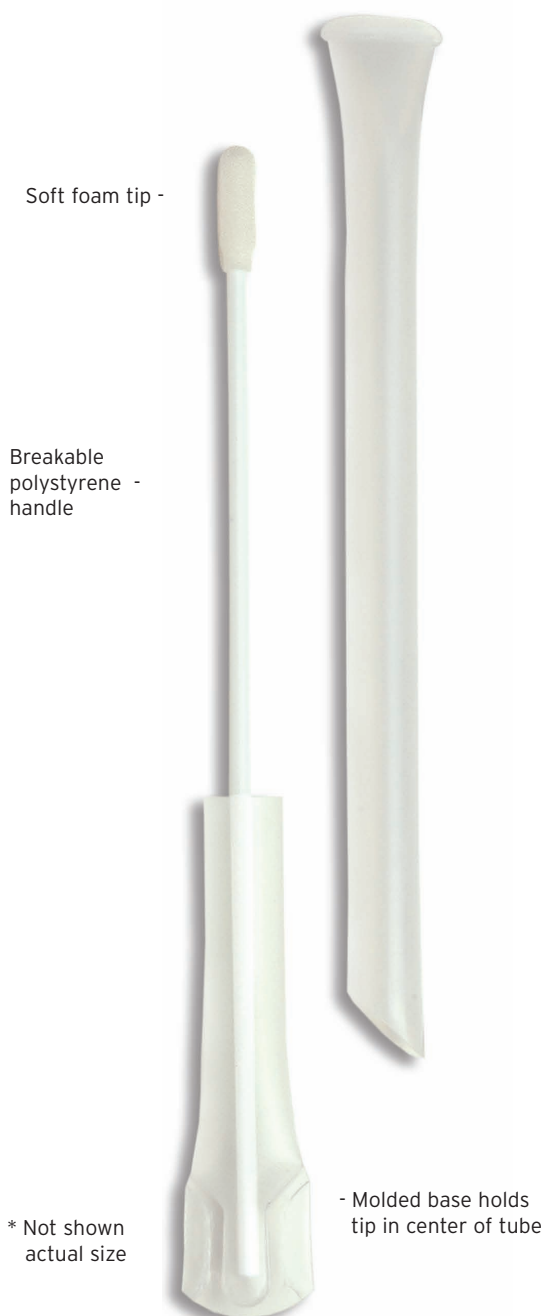

Description

This specialty applicator is tipped with 100ppi medical grade polyurethane foam. It is bonded securely to a hollow polystyrene handle that snaps when bent at a 90° angle. It comes packed in a dry transport device for ease in transport to the lab. The molded base holds tip in center of tube to help avoid loss of sample material during transport.

Dimensions

5.979" Overall Length (151.867mm)

Tip

■ 0.188" Dia. (4.775mm)

■ 0.63" Length (16mm)

Shaft

■ 0.100" Dia. (2.54mm)

■ 5.875" Length (149.23mm)

Standards

- CE marked for medical use
- ISO 9001:2008
- ISO 13485:2003
- FDA GMP's
- Certificate of sterility and/or compliance available with each shipment upon request

Qualities

- Medical grade quality
- Foam has excellent collection and release properties
- Can be used dry or with media
- Soft tip for patient comfort
- Packed in a reclosable dry transport tube for ease in transport to the lab

Applications

- Buccal swab
- For general specimen collection for testing and screening
- Used in Pediatrics, Ophthalmology, Veterinary care, and Oropharyngeal applications
- DNA testing

Ref #	Packaging
25-1506 1PF TT MC (sterile)	10/100/1 (10 boxes of 100 tubes per case, 1 per tube)



Puritan Medical Products Company LLC
Guilford, Maine 04443-0149 USA

Tel 800-321-2313 (US and Canada Only)
207-876-3311

Fax 800-323-4153
207-876-3130

Web www.puritanmedproducts.com

Email sales@puritanmedproducts.com

1-800-321-2313
puritanmedproducts.com